

Product Specifications

Physical Specifications

Dimensions	240 x 160 x 85 mm (L x W x H)
Weight	< 1.5 kg (approx.)
Display size	2.4 inch LCD screen

Performance specifications

Modes	CPAP, APAP
Pressure Range	4 - 20 cmH2O
Ramp Time	0 - 45 min / Auto
Flex	Off, Level 1 - 3
Automatic Therapy Start/Stop	Available
Humidifier Levels	Level 0 – 6
Tubing compatibility	15mm, 22mm
Filters	Standard, Ultra-fine (optional)
Altitude Compensation	Automatic
IP protection	IP22
Sound Pressure Level	<26 dB
External Communication	via Bluetooth

Power Supply

AC Power Consumption	100 - 240 V AC, 50 / 60 Hz, 1.0 A max
Device Input	24 V, 3.75 A

Warranty

2 years from the date of shipment



vREM

CPAP Device

Makes you sleep with comfort








vTitan Corporation Private Limited

New Woodlands Building, Estancia IT Park,
Plot No. 140,151 Vallancheri Village,
Chengalpattu District,
Tamil Nadu - 603202, India.

www.vtitan.com
sales@vtitan.com

©2024 vTitan Corporation Pvt. Ltd.
980-1065 Rev 04

Key Features

-  2.4" LCD with user friendly interface
-  Easy menu access
-  Check Mask fit feature to reduce leakage
-  Bluetooth application for instant Data Transfer
-  Visual and audible alarms

Comfort Features



Auto Adjusting Algorithm

Adjusts the therapy pressure automatically on a breath by breath basis.



Auto Ramp

Begins the therapy at a lower pressure and automatically increases once the user is asleep. Enhances patient comfort.



Flex Feature

Reduces the therapy pressure during exhalation and makes the user comfortable to breathe out.



Automatic Therapy Start/Stop


Automatically initiates the therapy once the device senses a normal breath and stops the therapy when the user removes the mask.



Integrated Humidifier

Adjusts the humidity levels during the therapy and helps to reduce dry nose and sore throat.

vREM Mobile Application

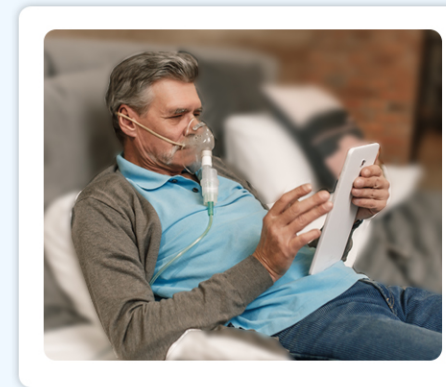

- View and track the therapy data after each session.
- Feedback about the treatment and the therapy results.
- Personalize the treatment with the help of alerts, reminders and personal goals.
- Review, store and share the reports externally.
- Visualize detailed graphical representation by sharing therapy data with OSCAR software