

Product Specifications

Physical Specifications

Dimensions 180 x 160 x 65 mm (L x B x H)

Weight 1 kg

Performance Specifications

Flow Rate 0.1 – 50 ml/s

Void Volume Up to 1000 ml

Volume Accuracy $\pm 2\%$ of voided volume

Modes Manual, Auto

Parameters Measured Maximum flow rate, Average flow rate, Voided volume, Voiding time, Flow time, Time to max, Hesitancy time

Sensor Load cell with 3 kg capacity

Alarms Jar displaced, Device overloaded, Device disconnected, Battery low, Battery depleted, Sensor failure

Ingress Protection IPX4

Data Communication Bluetooth

Battery

Type 4 AA batteries

Backup 10 days

Recharging Time 6 hrs

Accessories

Standard Charger, Uroflowmeter stand, Calibration stone

Optional Funnel, Collecting jar, Micturition chair

Warranty 1 year from the Date of Shipment

vTitan Corporation Private Limited

New Woodlands Building, Estancia IT Park,
Plot No. 140,151 Vallancheri Village,
Chengalpattu District,
Tamil Nadu - 603202, India.
Toll Free Number: 18005721664

www.vtitan.com

sales@vtitan.com

©2025 vTitan Corporation Pvt. Ltd.

980-1082 Rev 08











Accuflow® Uro


Uroflowmeter

- ◆ Compact & Portable
- ◆ Enhanced safety & simplicity
- ◆ Hassle-Free Operation
- ◆ Wireless Data Acquisition
- ◆ Real - time Display
- ◆ User-friendly

Accuflow Uro Key Features

-  Lightweight
-  High precision load cell sensor
-  Adjustable Urine Density and Flow Scale
-  Pre-configurable diagnostic codes
-  Wireless and Real-Time Data Acquisition
-  Sharing of report via email and other applications
-  Generation of detailed reports with patient data, results and nomograms
-  Multi-platform software for different use environments

User-Intuitive Application

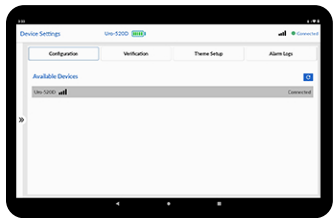

Automatic Device Configuration

Open the app to automatically pair with the device via Bluetooth




Effective Verification of Calibration

Provision for verifying the calibration to ensure high performance accuracy




Auto / Manual Modes of Operation

Auto Start / Stop feature that assures ease of use
Manual mode to calculate hesitancy

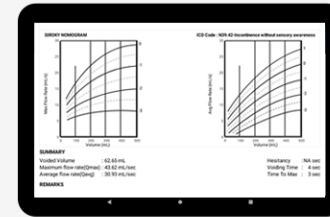


Comprehensive Reports

View detailed reports for current as well as previous tests
Wirelessly print or share the generated report via e-mail or other applications

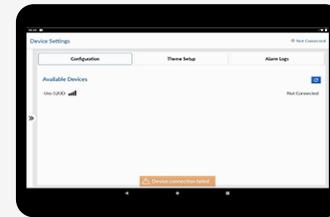

Prompt Visual Alarms

Lucid alarm indications on the app to alert the user at the right time

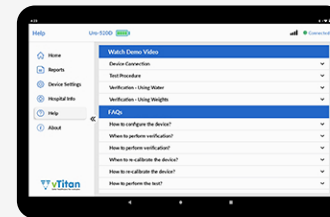

Help Section

Watch demo videos, view FAQs, Instructions For Use and other resources to know more about the device